

Portable Box-Type Automatic Wrapping Toilet

Wrappon·Green

Quick Reference Guide


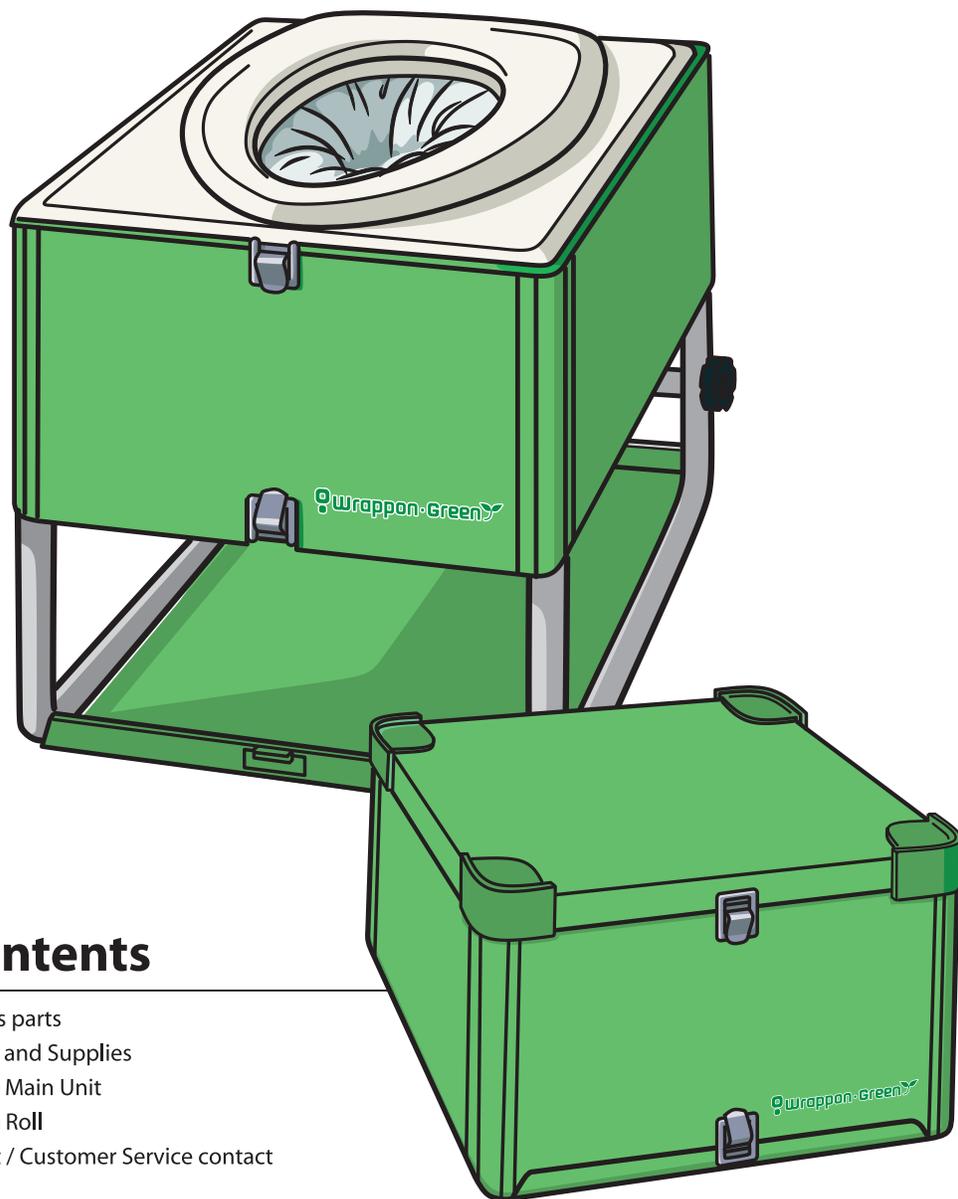


Table of Contents

- P.1 • • • Name of various parts
- P.3 • • • Accessories List and Supplies
- P.4 • • • Assembling the Main Unit
- P.5 • • • Installing a Film Roll
- P.7 • • • Using the Toilet / Customer Service contact


Name of various parts

【Front】

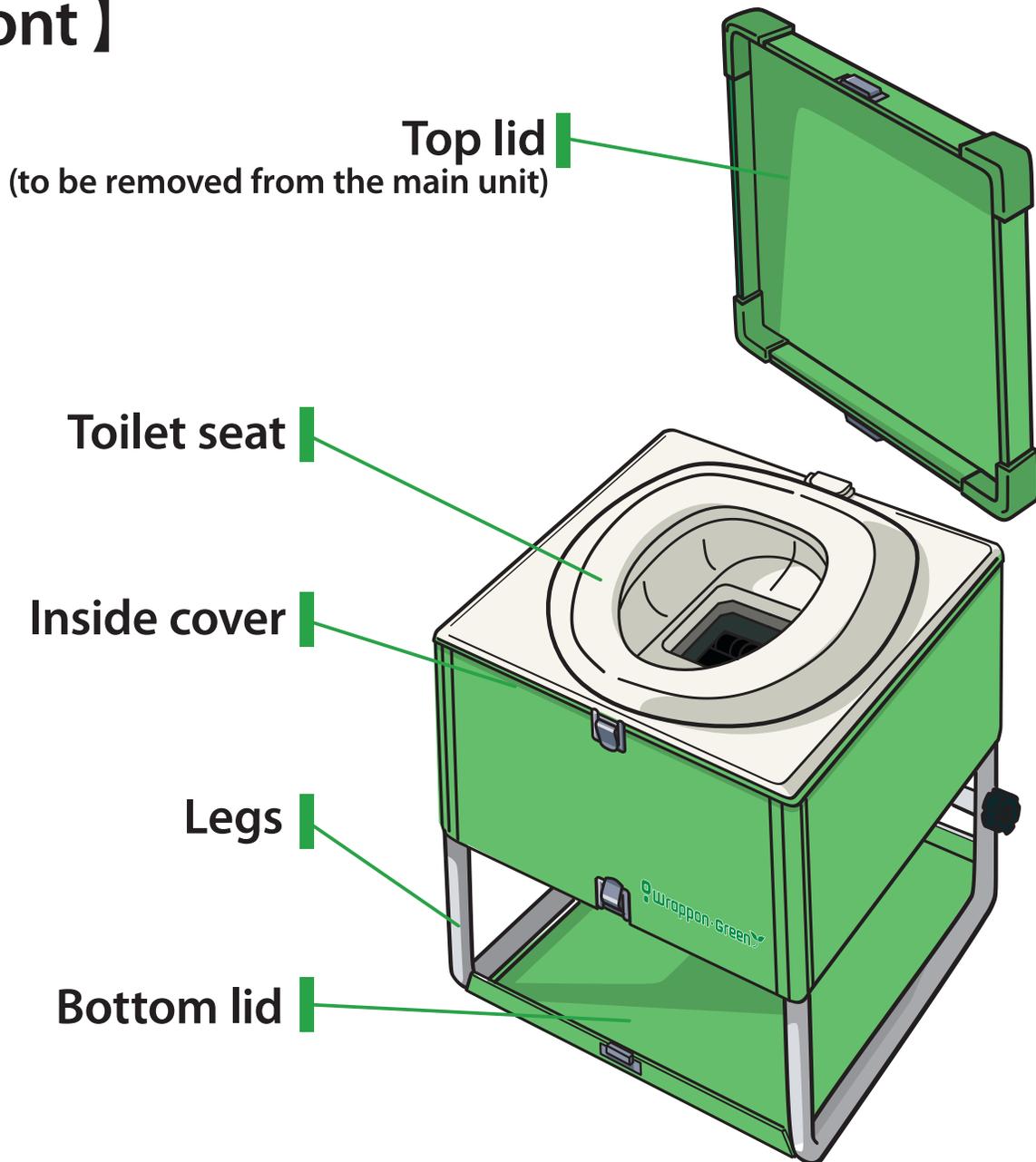

【AC adapter】

Power supply cord
plug-in port

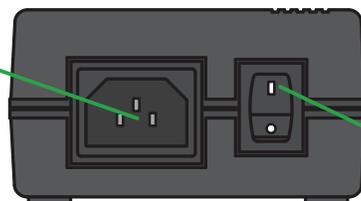


Main power
supply

【Rear】

Reinforcement bar

Screws

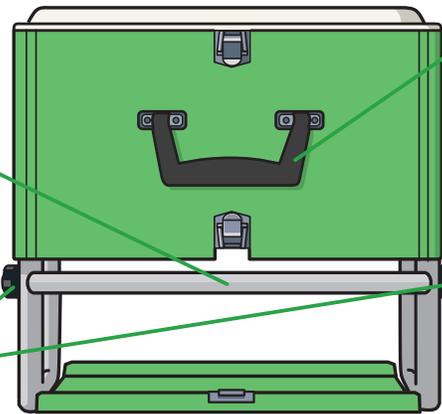


Handle

【Inside】 Top view with the toilet seat removed

Yellow switch
(initialization)

Blue switch
(for film feeding)

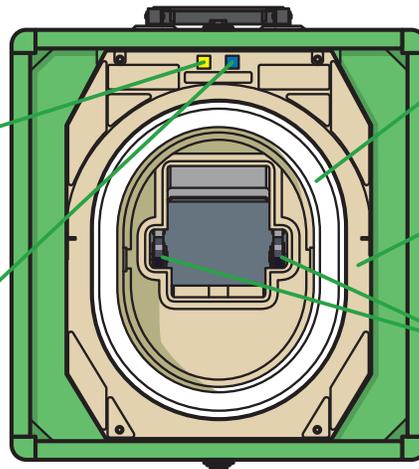

Film guide

Top plate


Roller

【Remote Control Unit】

(Hook-and-loop tape provided on back)

Display window

Button

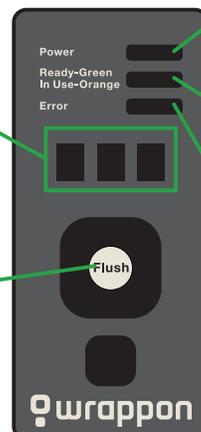


Power supply
indicator
(green)

Operation indicator
(green/orange)

Error indicator
(red)

Accessories List


Before using the product, check to make sure all of the accessories are intact.

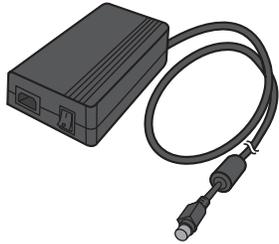

AC adapter
(Model No.: WRACFJ11S)



Power supply cord
(Model No.: WRPCFJ12U)

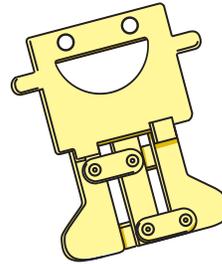

Film installation jig
(Model No.: COJ10001S)



User's Manual
(Model No.: GRMA0001U)

Quick Reference Guide
(Model No.: GRHA0001U)

Hook-and-loop tape
To be used for securing the Remote Control Unit.
Apply it to the desired location.

installed inside the product

Reinforcement bar
(Model No.: GRST0001S)

Two screws
Installed on the two ends of the reinforcement bar
(Model No.: GRBLN0B1S)

Inside cover (1)*
(Model No.: COMUUS02S)

Supplies


These supplies are required for using the product.
To order these, contact the seller where you purchased the product.



Paper coagulant
(Katamerusa, paper)
(For approximately 50 uses) 1 bag
(Model No.: COC0STP1U)














Polymer coagulant
(Katameru polymer)
(For approximately 50 uses) 1 box
(Model No.: COC0P001U)








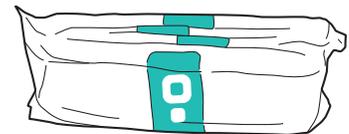





Film roll
(For approximately 50 uses) 1 roll
(Model No.: COF150C1U)

Assembling the Main Unit

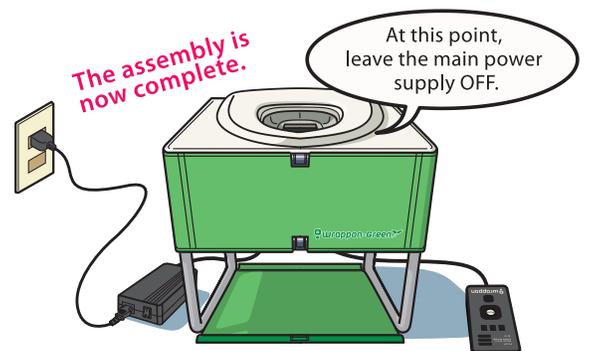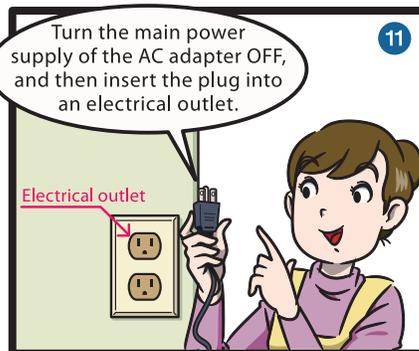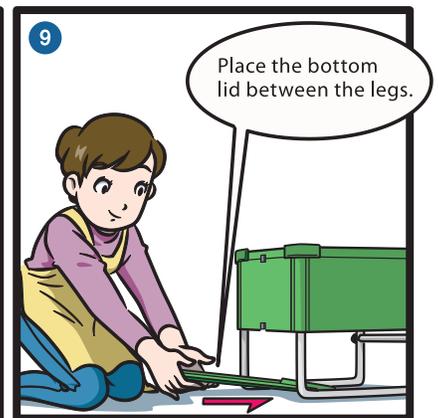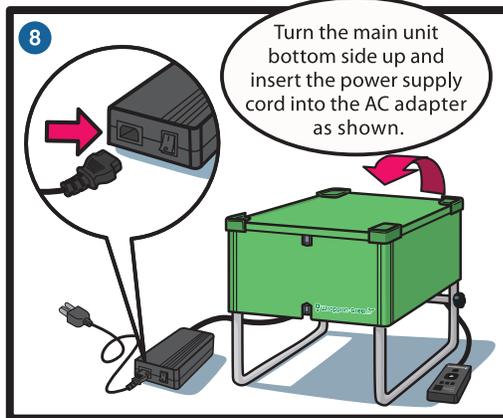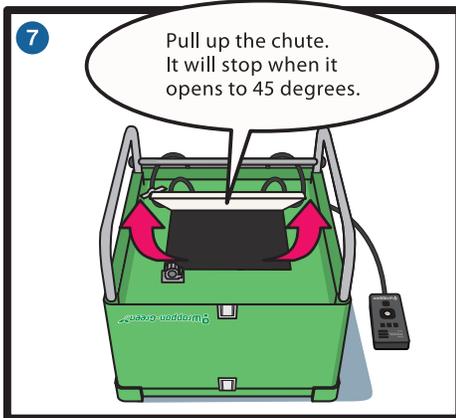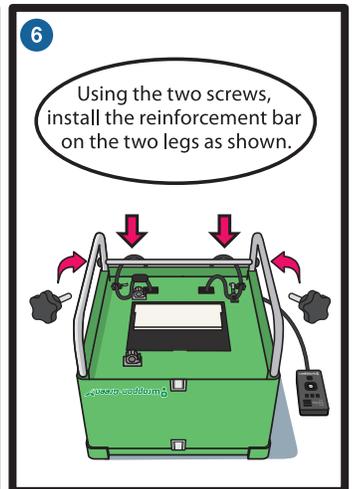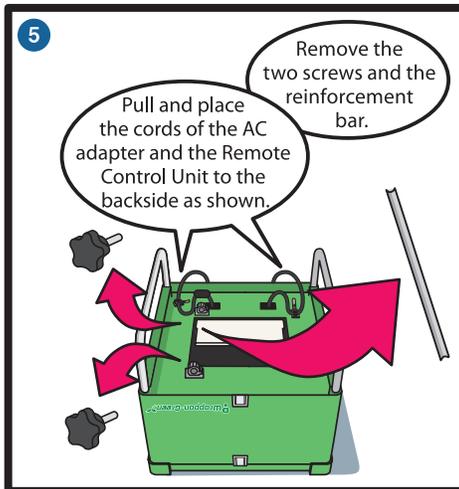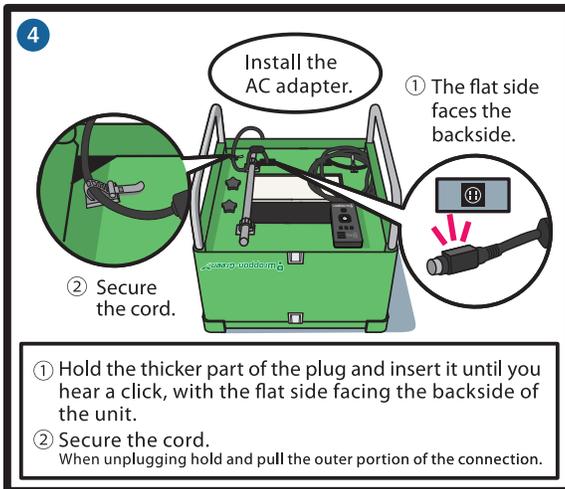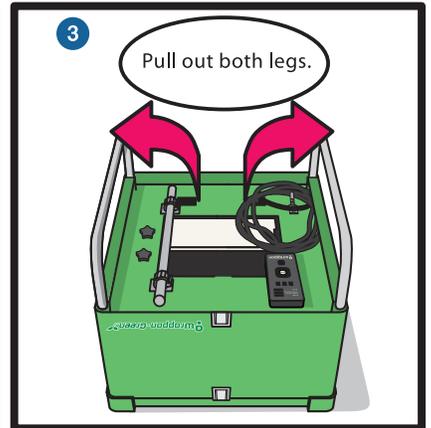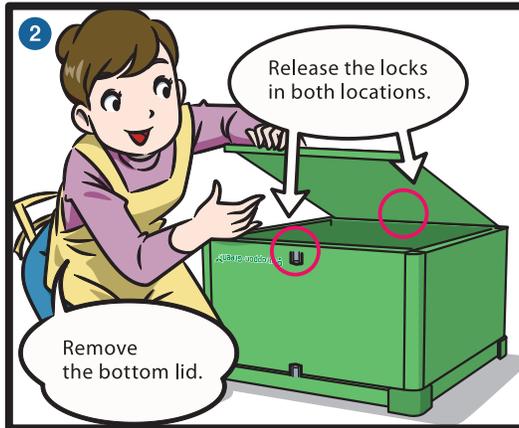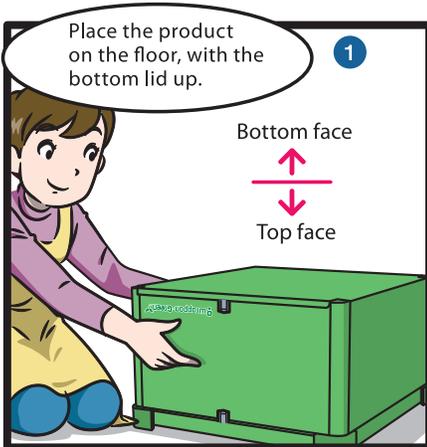

Installing a Film Roll

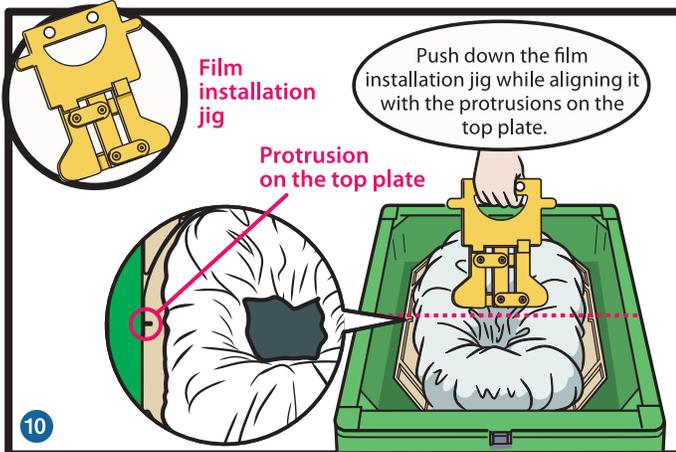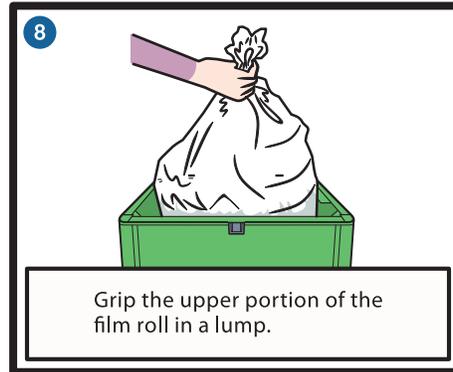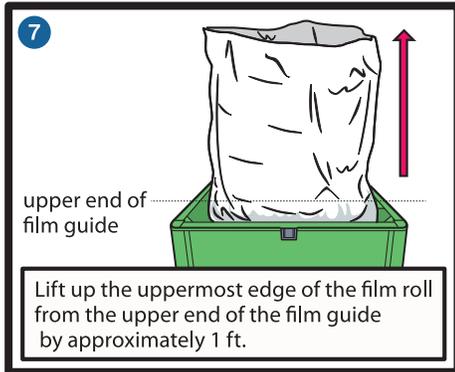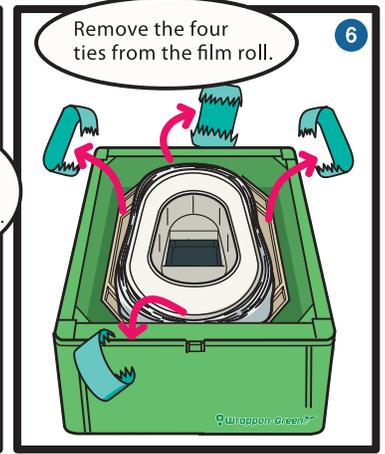
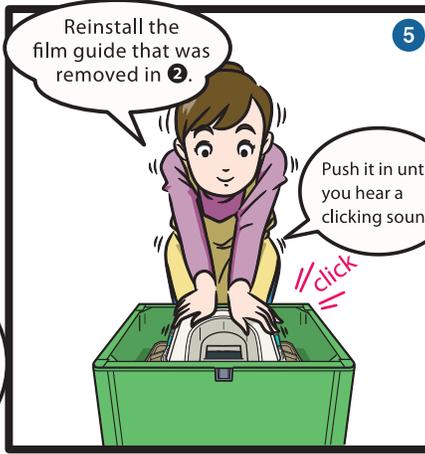
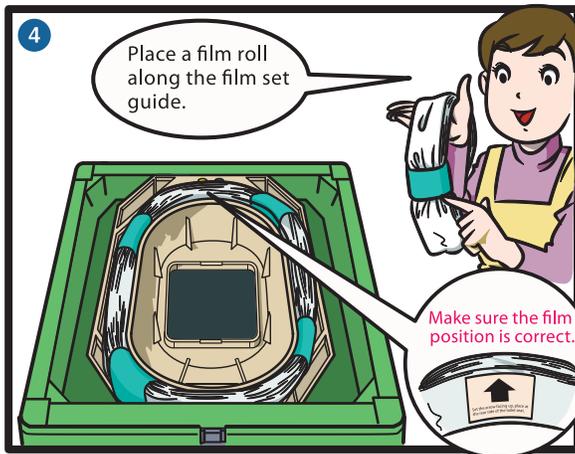
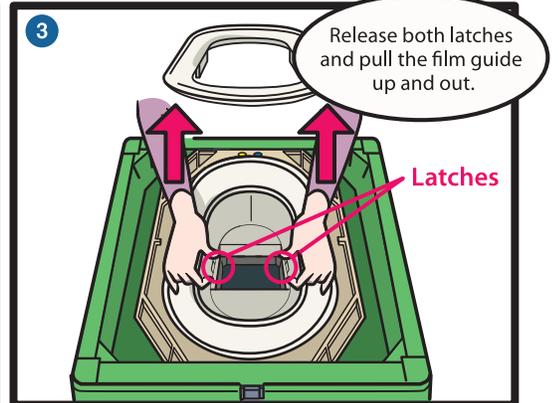
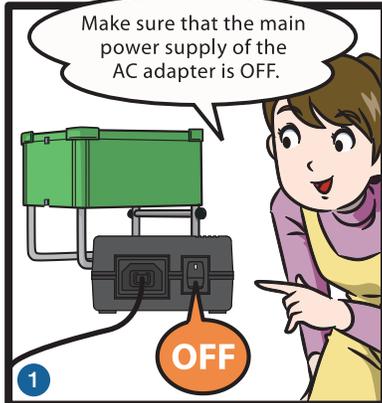


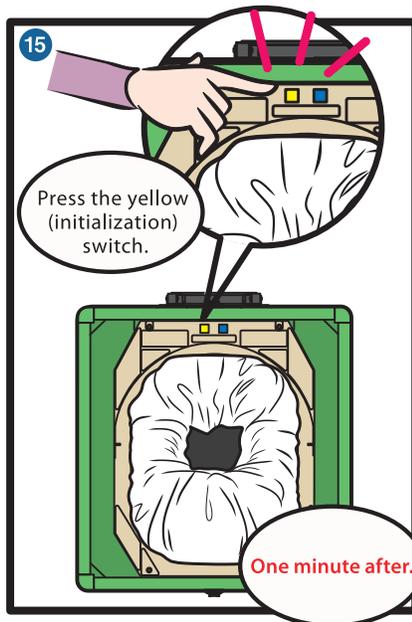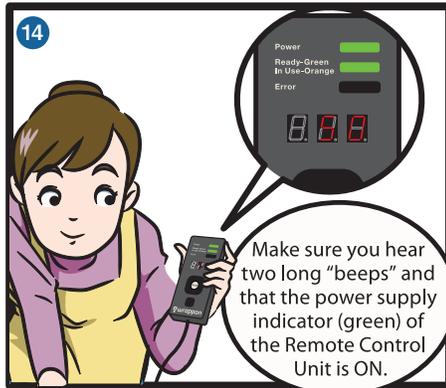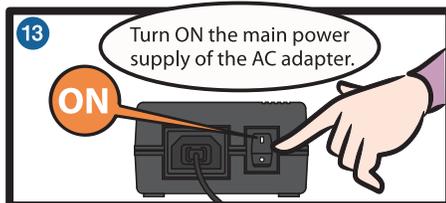

When the “STOP” mark appears on the film inside the toilet, it is time to replace the film roll.



When the “STOP” mark appears on the film

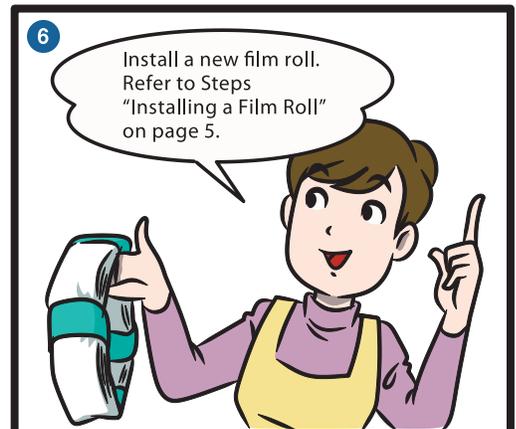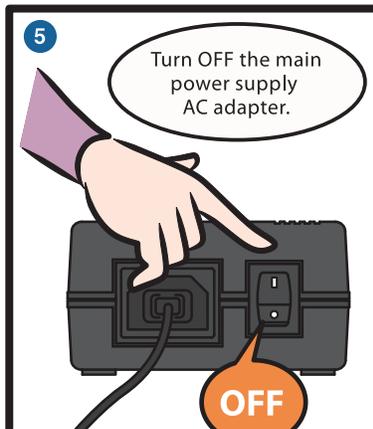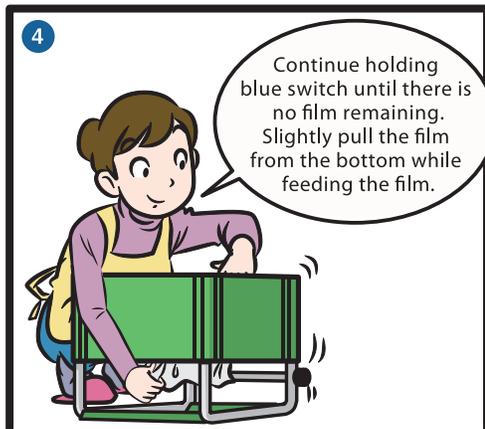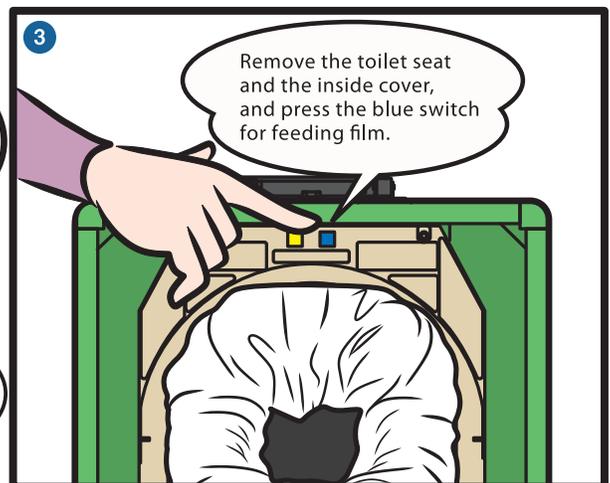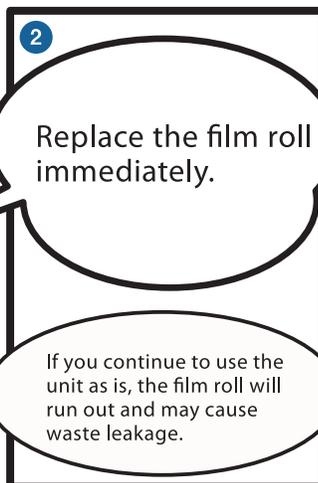

Using the Toilet

Do not put your hand in the toilet.
Sit on the toilet seat.

Checking before use.

Make sure that the operation indicator on the Remote Control Unit is lit green.
Check for the presence of a film roll.

